

# Backup

Wiki section

## Contents

- 1 Process to Identify Possible Acquisition Targets
  - ◆ 1.1 Advantage of Mergers and Acquisitions
  - ◆ 1.2 Methodology
  - ◆ 1.3 Company Profile
    - ◇ 1.3.1 Applied Medical Resources Corp.
    - ◇ 1.3.2 Bioteque Corp.
- 2 M&A Due Diligence Process
- 3 Phase 1: Landscape overview
  - ◆ 3.1 Ureteral Stent: Concept
  - ◆ 3.2 Background
  - ◆ 3.3 Ureteric stenting difficulties
    - ◇ 3.3.1 Encrustation
    - ◇ 3.3.2 Obstruction
    - ◇ 3.3.3 Infection
  - ◆ 3.4 The market need
  - ◆ 3.5 Intellectual property
    - ◇ 3.5.1 Search strategy
    - ◇ 3.5.2 Sample patents
- 4 Urinary Problems in men and women
- 5 Market Analysis
  - ◆ 5.1 Methodology
  - ◆ 5.2 Prevalence rate in US (women)
  - ◆ 5.3 Urinary incontinence severity among different age groups in US women
  - ◆ 5.4 Market potential for ureteral stent in US (women)
  - ◆ 5.5 Ureteral stent market forecast in US (women)
  - ◆ 5.6 Ureteral stent companies
    - ◇ 5.6.1 Ureteral stents of various companies
      - 5.6.1.1 Boston Scientific
      - 5.6.1.2 Cook Group
      - 5.6.1.3 Applied medicals
      - 5.6.1.4 Redi-Tech Medical Products
- 6 Like this report?
  - ◆ 6.1 Clinical Trials
    - ◇ 6.1.1 New trials
    - ◇ 6.1.2 Concluded trials
    - ◇ 6.1.3 Adverse Events
  - ◆ 6.2 Products
  - ◆ 6.3 Startup activity
- 7 Phase 2: Deeper Dive
  - ◆ 7.1 Scenario
  - ◆ 7.2 Deal analysis for a target company
  - ◆ 7.3 Design History File Review: Review components
  - ◆ 7.4 Sample report
    - ◇ 7.4.1 Performance/Functional Characteristics
  - ◆ 7.5 Potential DHF Review Outcomes
- 8 Phase 3: Post-acquisition integration
  - ◆ 8.1 Deadlines
  - ◆ 8.2 Documents and Ownership
- 9 Ureteral stents regulatory issues
  - ◆ 9.1 Pre-Market Notification
  - ◆ 9.2 Timeline Sheet
- 10 Like this report?
- 11 Contact Dolcera
- 12 Addition
  - ◆ 12.1 Patent Categorization: Interactive mind map linked to Dolcera Dashboard
  - ◆ 12.2 Product-Patent-Clinical Trials Mapping
  - ◆ 12.3 Product to Patent Mapping
  - ◆ 12.4 Dolcera Dashboard
- 13 Key Artifacts
- 14 Removed Sections

## Process to Identify Possible Acquisition Targets

Before identifying the possible acquisition targets, first companies needs to identify the objective of the merger and acquisition.

## Advantage of Mergers and Acquisitions

The most common motives and advantages of mergers and acquisitions are:-

- Accelerating a company's growth, particularly when its internal growth is constrained due to paucity of resources. Internal growth requires that a company should develop its operating facilities- manufacturing, research, marketing, etc. But, lack or inadequacy of resources and time needed for internal development may constrain a company's pace of growth. Hence, a company can acquire production facilities as well as other resources from outside through mergers and acquisitions. Specially, for entering in new products/markets, the company may lack technical skills and may require special marketing skills and a wide distribution network to access different segments of markets. The company can acquire existing company or companies with requisite infrastructure and skills and grow quickly.
- Enhancing profitability because a combination of two or more companies may result in more than average profitability due to cost reduction and efficient utilization of resources. This may happen because of:-
  - ◆ Economies of scale:- arise when increase in the volume of production leads to a reduction in the cost of production per unit. This is because, with merger, fixed costs are distributed over a large volume of production causing the unit cost of production to decline. Economies of scale may also arise from other indivisibilities such as production facilities, management functions and management resources and systems. This is because a given function, facility or resource is utilized for a large scale of operations by the combined firm.

- ♦ Operating economies:- arise because, a combination of two or more firms may result in cost reduction due to operating economies. In other words, a combined firm may avoid or reduce over-lapping functions and consolidate its management functions such as manufacturing, marketing, R&D and thus reduce operating costs. For example, a combined firm may eliminate duplicate channels of distribution, or create a centralized training center, or introduce an integrated planning and control system.
- ♦ Synergy:- implies a situation where the combined firm is more valuable than the sum of the individual combining firms. It refers to benefits other than those related to economies of scale. Operating economies are one form of synergy benefits. But apart from operating economies, synergy may also arise from enhanced managerial capabilities, creativity, innovativeness, R&D and market coverage capacity due to the complementarity of resources and skills and a widened horizon of opportunities.
- Diversifying the risks of the company, particularly when it acquires those businesses whose income streams are not correlated. Diversification implies growth through the combination of firms in unrelated businesses. It results in reduction of total risks through substantial reduction of cyclicity of operations. The combination of management and other systems strengthen the capacity of the combined firm to withstand the severity of the unforeseen economic factors which could otherwise endanger the survival of the individual companies.
- A merger may result in financial synergy and benefits for the firm in many ways:-
  - ♦ By eliminating financial constraints
  - ♦ By enhancing debt capacity. This is because a merger of two companies can bring stability of cash flows which in turn reduces the risk of insolvency and enhances the capacity of the new entity to service a larger amount of debt
  - ♦ By lowering the financial costs. This is because due to financial stability, the merged firm is able to borrow at a lower rate of interest.
- Limiting the severity of competition by increasing the company's market power. A merger can increase the market share of the merged firm. This improves the profitability of the firm due to economies of scale. The bargaining power of the firm vis-à-vis labour, suppliers and buyers is also enhanced. The merged firm can exploit technological breakthroughs against obsolescence and price wars.

**Note:** In this research report, researcher considered two following objective of merger and acquisition:

1. To Enter in emerging markets
2. To improve product portfolio

## Methodology

- **Step 1:** First, list all the players present in Ureteral Stent Market were identified. (The list of companies was retrieved from the FDA site from Registration & Listing database)
  - ♦ Number of companies that were identified in this step: 20
- **Step 2:** The second step involved identifying and eliminating companies that are large, and established players or the subsidiaries of the big players in the industry.
  - ♦ Number of companies eliminated in this step: 14
  - ♦ Number of (small) companies of interest left: 6
- **Step 3:** Once the large, and established players were eliminated, companies were compared based on various parameters and rated on the scale of 5 to identify the best target. Please check the following table and dashboard:
  - ♦ Number of potential target companies: 6

Company	Headquartered	Stent Details	Company Type	2010 Revenue (Mn)	No of Employees	Geographical Revenue Share	R&D investment (Mn)	No. of Patent	Technology Focus
Applied Medical Resources Corp.	USA	Mesh ureteral stent C-flex ureteral stent Silicone ureteral stent Tethered ureteral stent C-flex ureteral stent Ureteral stent Silhouette pediatric stent	Manufacturer	\$282	1900	USA: 62% Asia: 32% Rest of World: 6%	\$38.3 (13%)	212	High
Bioteque Corp.	Taiwan	Bioteque Double Pigta	Contract Manufacturer; Contract Sterilizer; Manufacturer	\$27.85	350	Taiwan: 100%	\$1.4 (5%)	16	Medium
Hobbs Medical, Inc.	USA	HM Ureteral Double Pigtail Stent	Manufacturer	\$2.30	22	USA: 97.2% Japan: 2.3% Others: 0.5%	\$0.06 (2.87%)	1	Low
Lake Region Medical Limited	Ireland	M-Wires	Contract Manufacturer	\$84.86	611	Europe: 68% USA: 21% Others: 11%	\$8.48 (3.85%)	8	Medium
	USA			\$23.00	180			0	Low

Martech Medical Products		Ureteral Stent	Contract Manufacturer			USA: 78.3% Europe: 12.7% Others: 9%	\$1.1 (4.79%)		
Allium Medical	Israel	URS - Ureteral Stent TPS - Triangular Prostatic Stent BUS - Bulbar Urethral Stent RPS - Round Posterior Urethral Stent BIS - Biliary Stent	Manufacturer	\$2.04	20	Domestic: 88.6% International: 11.4%	\$0.082 (4%)	2	Low

**Note:** Data for Geographical segment, R&D expenditure and Technology focus has been masked.

This is supposed to be a flash animation. You'll need the flash plugin and a browser that supports it to view it.

- **Step 4:** Finally, after an in-depth analysis of the potential target companies? on various parameters, following companies identified to be best possible target companies.
  - ♦ To improve Product Portfolio : **Applied Medical Resources Corp.**
  - ♦ To enter in Emerging Markets : **Bioteque Corp.**

## Company Profile

### Applied Medical Resources Corp.

APPLIED MEDICAL RESOURCES CORP.	
<b>Revenues</b>	2010: \$280 million
<b>Net profit (2010)</b>	\$20.89 million
<b>R&amp;D Investment</b>	2010: \$ 38.3 million (13% of revenues)
<b>Number of employees</b>	1900
<b>Year Established</b>	1987
<b>Headquarters</b>	USA
<b>Key People</b>	CEO: Said S. Hilal
<b>Products &amp; Technology</b>	Specialty Areas:Cardiac/Vascular, Colorectal, GYN, Urology New Products:Epix (Laparoscopic Instrumentation), GelPOINT (Advanced Access Platform), Kii Fios (First EntrySystem)
<b>Products in Ureteral Stent</b>	C-flex ureteral stent, Silicone ureteral stent, Tethered ureteral stent, C-flex ureteral stent, Ureteral stent Silhouette pediatric stent
<b>Geographical revenue breakdown (2010)</b>	USA: \$182.70 million (65%) Others: \$99.47 million (35%)
<b>Company Overview</b>	Applied Medical Resources Corporation is a new generation medical device company founded in 1987 and headquartered in Southern California. It is involved in developing, manufacturing and marketing of innovative products for Minimally Invasive Surgery, Cardiac, Vascular, Urological, Colorectal, Bariatric, Obstetric, Gynecologic and General Surgery. The product portfolio covers 25 technologies and more than 700 products. The company has spread its business globally across 75 countries including Africa, Middle East,Americas, Caribbean, Asia, Australia and Europe through its network of international distributors.

Bioteque Corp.	
<b>Revenues</b>	2010: \$27.85 million
<b>Net profit (2010)</b>	\$20.89 million
<b>R&amp;D Investment</b>	2010: \$38.3 million (13% of revenues)
<b>Number of employees</b>	350
<b>Year Established</b>	1991
<b>Headquarters</b>	Taiwan
<b>Key People</b>	
<b>Products &amp; Technology</b>	NEPHROLOGY, UROLOGY, RADIOLOGY, CARDIOLOGY, RESPIRATORY CARE CRITICAL CARE, IV ADMINISTRATION THERAPY, MOLDING PARTS
<b>Products in Ureteral Stent</b>	C-flex ureteral stent, Silicone ureteral stent, Tethered ureteral stent, C-flex ureteral stent, Ureteral stent Silhouette pediatric stent
<b>Geographical revenue breakdown (2010)</b>	NA
<b>Company Overview</b>	Bioteque Corporation manufactures and sells medical devices in Taiwan. It offers medical disposables for use in hemodialysis access, endovascular treatment, and other fields. It offers blood tubing lines, AVF needles, transducer protectors, and on line HDF with check valves or without check valves; IV infusion bags, precision IV infusion sets, drainage bags, insufflation tubing sets/filters, and various surgical drainage tubes; a range of medical components, which comprise blood tubing line components, percutaneous drainage components, infusion bag components, AVF needle components, precision IV infusion set components, and IV infusion bag components; a range of thermoplastic polyurethane catheters, including pigtail drainage catheter sets, double pigtail ureteral stent sets, biliary drainage catheters, percutaneous nephrostomy kits, and dialysis catheters; and other medical disposable products, such as closed suction catheters and artificial nose. The company also provides endovascular products, which consist of percutaneous transluminal coronary angioplasty, percutaneous transluminal angioplasty, angiography catheters, guiding catheters, sheath introducers, MRI/CT/angiography syringes, micro catheters, and hydrophilic coated guidewires. Bioteque Corporation was founded in 1991 and is based in Taipei, Taiwan.

## M&A Due Diligence Process

### Phase 1: Landscape overview

#### Ureteral Stent: Concept

An antimicrobial ureteral stent, which inhibits encrustation and bacterial colonization while maintaining patient comfort.

- Ureteral stent: resists migration, resists fragmentation, is kink resistant and radiopaque.
- Bacterial colonization: antimicrobial activity for up to two weeks.
- Patient Comfort: stent has a low coefficient of friction (value) for ease of insertion and will soften on implant at body temperature to maintain patient comfort.

#### Background

Ureteral stents are used in urological surgery to maintain patency of the ureter to allow urine drainage from the renal pelvis to the bladder. These devices can be placed by a number of different endourological techniques. They are typically inserted through a cystoscope and may also be inserted intraoperatively. Indwelling ureteral stents help to reduce complications and morbidity subsequent to urological and surgical procedures. Frequently, ureteral stents are used to facilitate drainage in conjunction with Extracorporeal Shock Wave Lithotripsy (ESWL) and after endoscopic procedures. They are also used to internally support anastomoses and prevent urine leakage after surgery. Ureteral stenting may almost eliminate the urological complications of renal transplantation.


The advent of ESWL and the more recent barrage of endourological techniques have increased the indications for ureteral stents (Candela and Bellman 1997). Indications for use include:

- Treatment of ureteral or kidney stones
- Ureteral trauma or stricture
- Genitourinary reconstructive surgery
- Hydronephrosis during pregnancy

- Obstruction due to malignancy
- Retroperitoneal fibrosis

The need for ureteral stents range from a few days to several months. For patients with serious urological problems, ureteral stent maintenance may become a life-long necessity. Unfortunately, there are many problems associated with using ureteral stents.

## Ureteric stenting difficulties

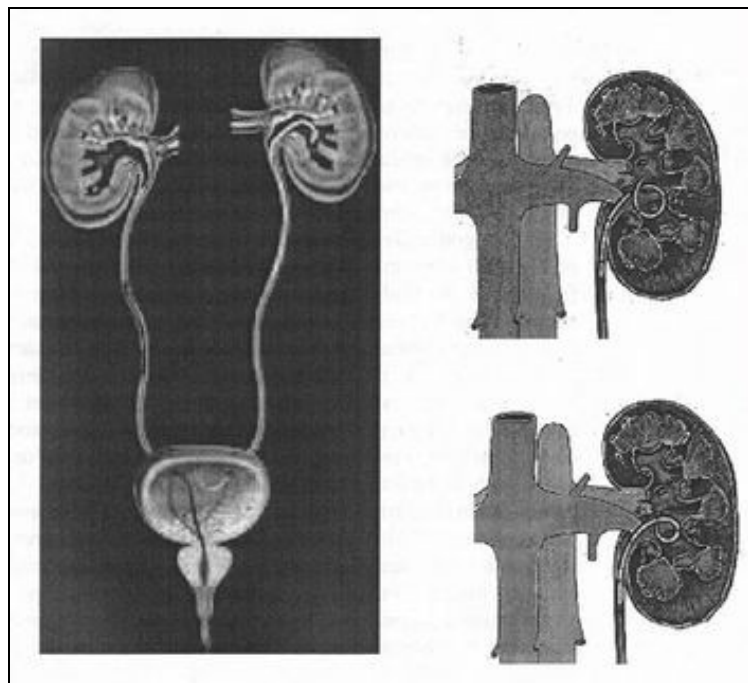

Double-J and Pigtail ureteral stents

Common	Rare
<ul style="list-style-type: none"> <li>• Trigonal irritation</li> <li>• Haematuria</li> <li>• Fever</li> <li>• Infection</li> <li>• Tissue inflammation</li> <li>• Encrustation</li> <li>• Biofilm formation</li> </ul>	<ul style="list-style-type: none"> <li>• Obstruction</li> <li>• Kinking</li> <li>• Ureteric rupture</li> <li>• Ureteric perforation</li> <li>• Stent misplacement</li> <li>• Stent migration</li> <li>• Stent misfit</li> <li>• Stent forgotten</li> <li>• Tissue hyperplasia</li> </ul>

Today, elastomeric materials, such as silicones, polyurethanes and hydrogel-coated polyolefins are used, with no clear winner, which can withstand the urinary environment.

- Although silicone has better long-term stability than other stent materials, its extreme flexibility makes it difficult to pass over guidewires and through narrow or tortuous ureters.
- Polyethylene is stiffer and easier to use for patients with strictures; however, it has been known to become brittle with time leading to breakage and is no longer commercially available. \* Polyurethane has properties that fall in between polyethylene and silicone; however, stent fracture also has been an issue with polyurethanes.

Attempts have been made to develop polymers with a combination of the best of all properties. The key players are C-Flex (Concept Polymer Technologies), Silitek and Percuflex (Boston Scientific).

- C-Flex is proprietary silicone oil and mineral oil interpenetrated into a styrenelolefin block copolymer with the hope of reduced encrustation.
- Silitek (Medical Engineering Corporation) is another silicone-based copolymer.
- Percuflex is a proprietary olefinic block copolymer.

Metallic stents have been used recently to treat extrinsic ureteric obstructions. The effect of synthetic polymers on the urothelium of the urinary tract seems to be dependent on the bulk chemical composition of the polymer, the chemical composition of its surface, coatings on the device surface, smoothness of the surface and coefficient of friction.

Typically, most ureteral stents are made of relatively smooth catheters. Koleski et al., (2000) tested a longitudinally grooved ureteral stent made by Circon in the pig ureter. The results indicated that the grooved stent led to better drainage than a conventional stent. Their opinion is that the ureter wall has a better chance of collapsing over a smooth surface than a grooved surface, especially when debris is present. Stoller (2000) had the same experience with the SpiraStent(Urosurge Corp.). This helical stent was superior at passing stones than a conventional smooth stent.

There are a variety of ureteral stent configurations with different anchoring systems. Most stents today have a double pigtail anchoring system. (Tolley, 2000), Dunn et al. (2000) conducted a randomized, single-blind study comparing a Tail stent (proximal pigtail with a shaft which tapers to a lumenless straight tail) to a double pigtail stent. The Tail stent was found to be better tolerated than the double-pigtail concerning lower urinary tract irritative symptoms. A double-J ureteral stent and a flexible ureteropyeloscope are shown in the first diagram. The other two diagrams show a pigtail ureteral stent in place; the end of the pigtail is facing away from the ureteral opening in the second of these two diagrams.

Early adverse effects of ureteral stenting include lower abdominal pain, dysuria, fever, urinary frequency, nocturia and hematuria. Patient discomfort and microscopic hematuria happen often. Major late complications include stent migration, stent fragmentation or more serious hydronephrosis with flank pain and infections.

Late complications occurred in one third of the patients in a prospective study using both silicone and polyurethane double pigtail stents (110 stents) in 90 patients. Stent removal was necessary in these patients. Others also have found this percentage of late complications. Device-related urinary tract infection and encrustation can lead to significant morbidity and even death and are the primary factors limiting long-term use of indwelling devices in the urinary tract. Microbial biofilm and encrustation may lead to stone formation. This is typically not a problem when stents are used for short-term indications. Problems of biofilm formation, encrustation and stent fracture occur in patients with long-term indwelling stents.

Typically, manufacturers advise periodic stent evaluation. Cook polyurethane stent removal is recommended at 6 months and 12 months for silicone (Cook product literature). However, stents that are intended for long-term use are usually changed at regular intervals, as frequently as every 3 months.

Forgotten stents are a problem. Monga et al., 1995 found that 68% of stents forgotten more than 6 months were calcified and 10% were fragmented. Multiple urologic procedures were necessary to remove the stones. Long-term effects of these forgotten stents may lead to voiding dysfunction and renal insufficiency. Schlick, et al., 1998 are developing a biodegradable stent that will preclude the need for stent removal.

### **Encrustation**

The urinary system presents a challenge because of its chemically unstable environment. Long-term biocompatibility and biodegradability of devices have been problems due to the supersaturation of uromucoids and crystalloids at the interface between urine and the device. Encrustation of ureteral stents is a well-known problem, which can be treated easily if recognized early. However, severe encrustation leads to renal failure and is difficult to manage (Mohan-Pillai et al., 1999). All biomaterials currently used become encrusted to some extent when exposed to urine.

The encrusted deposits can harbor bacterial biofilms. In addition, they can render the biomaterial brittle which causes fracture in-situ, a serious problem especially associated with the use of polyethylene and polyurethane ureteral stents (although silicone stents have also been reported to fracture). Stent fragments can migrate to the bladder or renal pelvis with serious repercussions.

Surface science techniques were used to study three stent types after use in patients. The stent type, duration of insertion and age or sex of the patient did not correlate significantly with the amount of encrustation (Wollin et al., 1998). However, it has been suggested that factors which affect the amount of encrustation include the composition of the urine, the type of invading and colonizing bacteria and the structure and surface properties of the biomaterial used (Gorman 1995). A low surface energy surface seems to resist encrustation compared with a high surface energy surface (Denstedt et al., 1998).

Many different types of stone can form in the urinary tract. Calcium oxalate, calcium phosphate, uric acid and cystine stones are metabolic stones because they form as a result of metabolic dysfunction. They usually are excreted from the urinary tract. Struvite (magnesium ammonium phosphate) and hydroxyapatite (calcium phosphate) are associated with infection (infection stones). These account for 1520% of urinary calculi. ESWL is used to break up the larger infection stones because they don't pass; recurrence of the problem occurs with incomplete removal. Infection stones can manifest as poorly mineralized matrix stones, highly mineralized staghorn calculi or as bladder stones which often form in the presence of ureteral stents. Urea-splitting bacteria colonize the surface and cause alkalization of the urine, which lowers the solubility of struvite and hydroxyapatite, and they deposit on the surface. Bacterial biofilm associated with encrustation is a common clinical occurrence. (Gorman and Tunney, 1997). It has been suggested that prevention of bacterial colonization would prevent encrustation because of their ultimate responsibility for its formation (Bibby et al., 1995).

An in vitro model was developed that produces encrustation similar to those seen in vivo (Tunney et al., 1996a). An experiment was conducted to compare the encrustation potential of various ureteral stent materials. The long-term struvite and hydroxyapatite encrustation of silicone, polyurethane, hydrogel-coated polyurethane, Silitek and Percuflex were compared. All of the materials developed encrustation, however, it was found by image analysis that the rates of encrustation varied on the different materials. Silicone had less encrustation (69% at 10 weeks) compared to the other materials (100%) at the same time point (Tunney et al., 1996b). Continuous flow models have also been developed which are more representative of conditions in the upper urinary tract. They are discussed by Gorman and Tunney, (1997). Efforts to reduce encrustation using new materials, smoother surfaces and hydrogel coatings have been attempted.

A hydrogel-coated C-flex stent (Hydroplus, Boston Scientific) was shown to have less epithelial cell damage and encrustation than other biomaterials and was recommended by the investigators for long-term use (Cormio, 1995). In addition, a poly(ethylene oxide)/polyurethane composite hydrogel (Aquavene, J & J) resisted intraluminal blockage in a urine flow model compared with silicone and polyurethane (Gorman et al., 1997a). Another advantage with Aquavene is that it is rigid in the dry state, which facilitates insertion past obstructions in the ureter and becomes soft on hydration providing comfort (Gorman and Tunney, 1997). Gorman et al. (1997b) concluded that the chance of stent fracture would be reduced if the ureteral stent side holes were eliminated. Urinary tract infection is another common major problem with the usage of ureteral stents. Initially, a conditioning film is deposited on the ureteral stent surface. The film is made up of proteins, electrolyte materials and other unidentified materials that obscure the surface properties of the stent material. Electrostatic interactions, the ionic strength and pH of the urine and differences in fluid surface tensions affect bacterial adhesion to the conditioning film. Subsequently, a microbial biofilm forms over time. The biofilm is composed of bacterial cells embedded in a hydrated, predominantly anionic mixture of bacterial exopolysaccharides and trapped host extracellular macromolecules.

### **Obstruction**

Obstruction of urine flow and urinary tract sepsis can result in continued growth of the biofilm. Colonization of devices implanted in the urinary tract can lead to dysfunction, tissue intolerance, pain, subclinical or overt infection and even urosepsis. Device related infections are difficult to treat and device removal is usually necessary. The biofilm has been found to impede the diffusion of antibiotics; in addition, the bacteria in the biofilm have a decreased metabolic rate, which also protects them against the effects of antibiotics (Wollin et al., 1998). Riedl, et al. (1999) found 100% ureteral stent colonization rates in permanent and 69.3% in temporary stents. Antibiotic prophylaxis did not prevent bacterial colonization and it was recommended that it not be used. On the other hand, Tieszer, et al. (1998) believe that fluoroquinolones can prevent infection. They also have found that some stents have denser encrustation than others, however, the stent material did not change the elements of the "conditioning film" adsorbed or alter its receptivity to bacterial biofilms.

### **Infection**

The predictive value of urine cultures in the assessment of stent colonization was examined in 65 patients with indwelling ureteral stents. It was found that a sterile urine culture did not rule out the stent itself being colonized (Lifshitz, et al., 1999). Patients with sterile urine culture may benefit from prophylactic antibiotics; however, the authors contended that the antibiotics must work against gram-negative uropathogens and gram-positive bacteria including enterococci. It is obvious that there is controversy in the literature whether prophylactic systemic antibiotics are useful with ureteral stent implant. However, antibiotics do not seem to prevent stent colonization. Denstedt et al. (1998) have found that ciprofloxacin, with a 3 day burst every 2 weeks, actually is adsorbed onto the stent which makes longer term treatment possible with reduced risk of bacterial resistance. There has been research targeted at coating or impregnating urinary catheters with antimicrobials and products are on the market, however, there are no antimicrobial ureteral stents approved by the FDA.

### **The market need**

It is clear that there is a need for a new material that will be able to resist encrustation and infection in the urinary tract. According to Merrill Lynch, ureteral stents represent an \$80 MM US market. Boston Scientific is in the lead with ~50% of the market followed by Maxxim (Circon), Cook and Bard is a smaller player. There are a number of other small contenders.

The use of ureteral stents is increasing; the indications for ureteral stenting have broadened from temporary or permanent relief of ureteric obstruction to include temporary urinary diversion following surgical procedures such as endopyelotomy and ureteroscopy and facilitation of stone clearance after ESWL (Tolley, 2000).



The use of ureteral stents for patients having ESWL for renal calculi is however controversial and seems to be related to the size of the stones and invasiveness of the procedure. According to survey results reported by Hollowell, et al. (2000), there is a significant difference in opinion concerning the use of stents with ESWL.

The number of ureteral stents used in patients with stones 2 cm or less treated with ESWL is significant in spite of the lack scientific evidence in support of this practice. Of 1,029 urologists returning surveys, for patients with renal pelvic stones 10, 15 or 20 mm treated with ESWL, routine stent placement was preferred by 25.3%, 57.1 % and 87.1 %, respectively. Urologists recommend using ureteroscopy rather than ESWL for distal ureteral calculi 5-10 mm.

## Intellectual property

### Search strategy

- Databases searched: US-G, US-A, EP-A, EP-B, WO, JP, DE, GB, FR
- Search scope: Title, Abstract or Claims
- Years: 1981-July 2008
- Search query: (ureter\* OR urether\* OR ureth\* OR uretr\*) AND (stent\*) AND (\*microb\* OR \*bacter\*)
- Results: **177 patents (82 unique patent families)**

### Sample patents

Patent	Assignee	Title	Abstract
<b>US6468649 B1</b>	SCIMED LIFE SYSTEMS INC	Antimicrobial adhesion surface	The present invention provides an implantable medical device having a substrate with a hydrophilic coating composition to limit in vivo colonization of bacteria and fungi. The hydrophilic coating composition includes a hydrophilic polymer with a molecular weight in the range from about 100,000 to about 15 million selected from copolymers acrylic acid, methacrylic acid, isocrotonic acid and combinations thereof.
<b>US5554147 A</b>	CApHCO, Inc.	Compositions and devices for controlled release of active ingredients	A method for the controlled release of a biologically active agent wherein the agent is released from a hydrophobic, pH-sensitive polymer matrix is disclosed and claimed. The polymer matrix swells when the environment reaches pH 8.5, releasing the active agent. A polymer of hydrophobic and weakly acidic comonomers is disclosed for use in the controlled release system. Further disclosed is a specific embodiment in which the controlled release system may be used. The pH-sensitive polymer is coated onto a latex catheter used in ureteral catheterization. A common problem with catheterized patients is the infection of the urinary tract with urease-producing bacteria. In addition to the irritation caused by the presence of the bacteria, urease produced by these bacteria degrade urea in the urine, forming carbon dioxide and ammonia. The ammonia causes an increase in the pH of the urine. Minerals in the urine begin to precipitate at this high pH, forming encrustations which complicate the functioning of the catheter. A ureteral catheter coated with a pH-sensitive polymer having an antibiotic or urease inhibitor trapped within its matrix will release the active agent when exposed to the high pH urine as the polymer gel swells. Such release can be made slow enough so that the drug remains at significant levels for a clinically useful period of time.
<b>US20030153983 A1</b>	SCIMED LIFE SYSTEMS INC	Implantable or insertable medical device resistant to microbial growth and biofilm formation	Disclosed are implantable or insertable medical devices that provide resistance to microbial growth on and in the environment of the device and resistance to microbial adhesion and biofilm formation on the device. In particular, the invention discloses implantable or insertable medical devices that comprise at least one biocompatible matrix polymer region, an antimicrobial agent for providing resistance to microbial growth and a microbial adhesion/biofilm synthesis inhibitor for inhibiting the attachment of microbes and the synthesis and accumulation of biofilm on the surface of the medical device. Also disclosed are methods of manufacturing such devices under conditions that substantially prevent preferential partitioning of any of said bioactive agents to a surface of the biocompatible matrix polymer and substantially prevent chemical modification of said bioactive agents

## Urinary Problems in men and women


- Both men and women have an increased risk for urinary incontinence as they get older, with men's rates rising steadily and women's rates peaking during menopause.
- The prevalence of incontinence in men of all ages is certainly lower than that for women.
- Women over 70, however, are twice as likely to have urinary incontinence as men of the same age.

Source: Urinary prevalence men Vs women

## Market Analysis

- We determined market data to have an idea about the market potential for ureteral stents.
- We have done this modeling for female population in US as women has the higher prevalence rate for urinary incontinence than men in all age groups.
- Prevalence increased with age, from 28% for 30- to 39-year-old women to 55% for 80- to 90-year-old women.
- 18% of respondents reported severe UI.
- The prevalence of severe UI also increased notably with age, from 8% for 30- to 39-year-old women to 33% for 80- to 90-year-old women.
- Among all, 9% reported slight UI, 15% reported moderate UI, 18% reported severe UI, and 58% reported no UI.

## Methodology

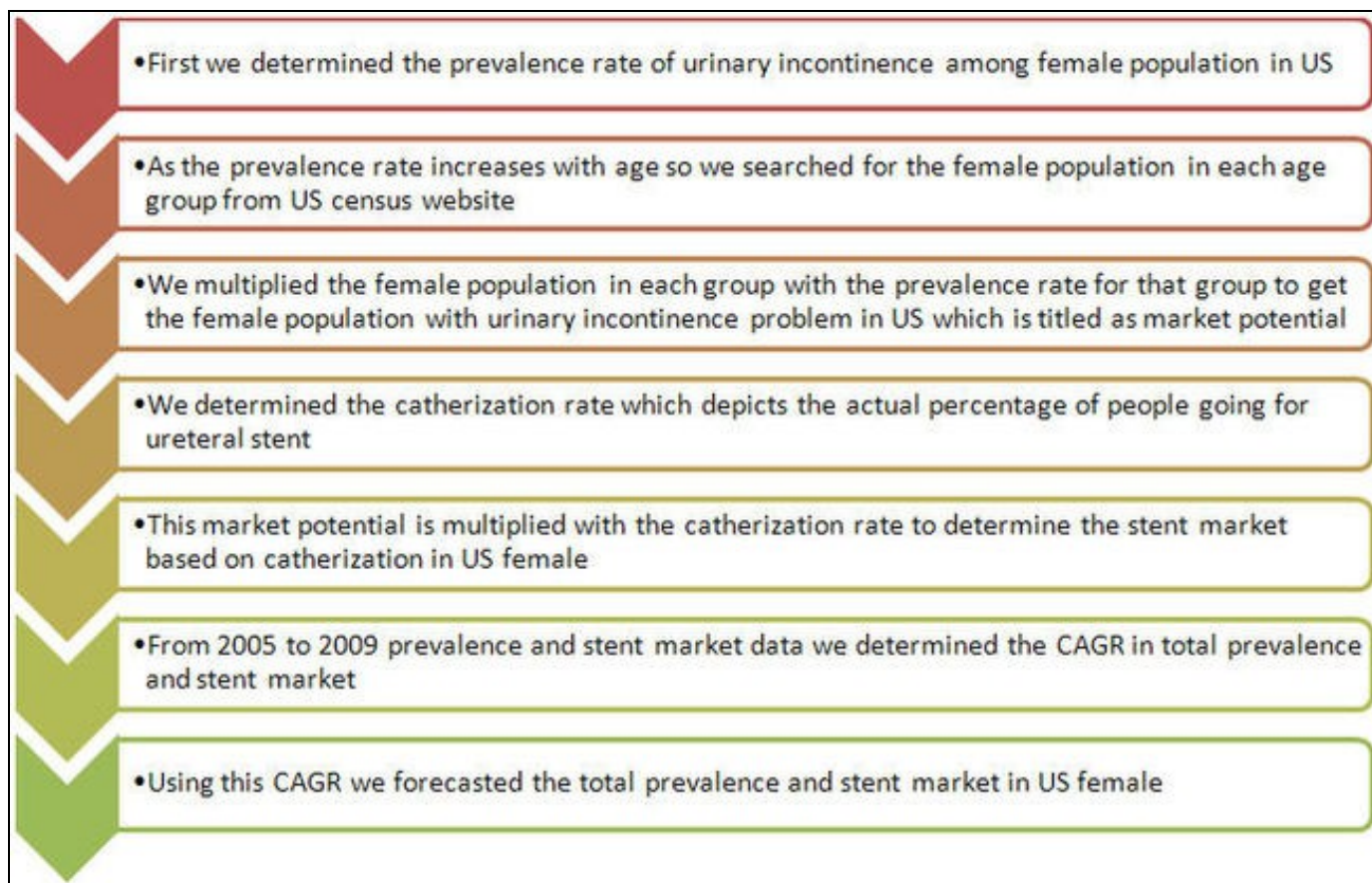


#### Prevalence rate in US (women)

Prevalence of Urinary Incontinence in US (women)	
Age (in yrs)	Population with Urinary incontinence (in %)
30-39	28%
40-49	41%
50-59	48%
60-69	51%
70-79	55%
80-90	54%

Source: Archives of internal medicine

#### Urinary incontinence severity among different age groups in US women



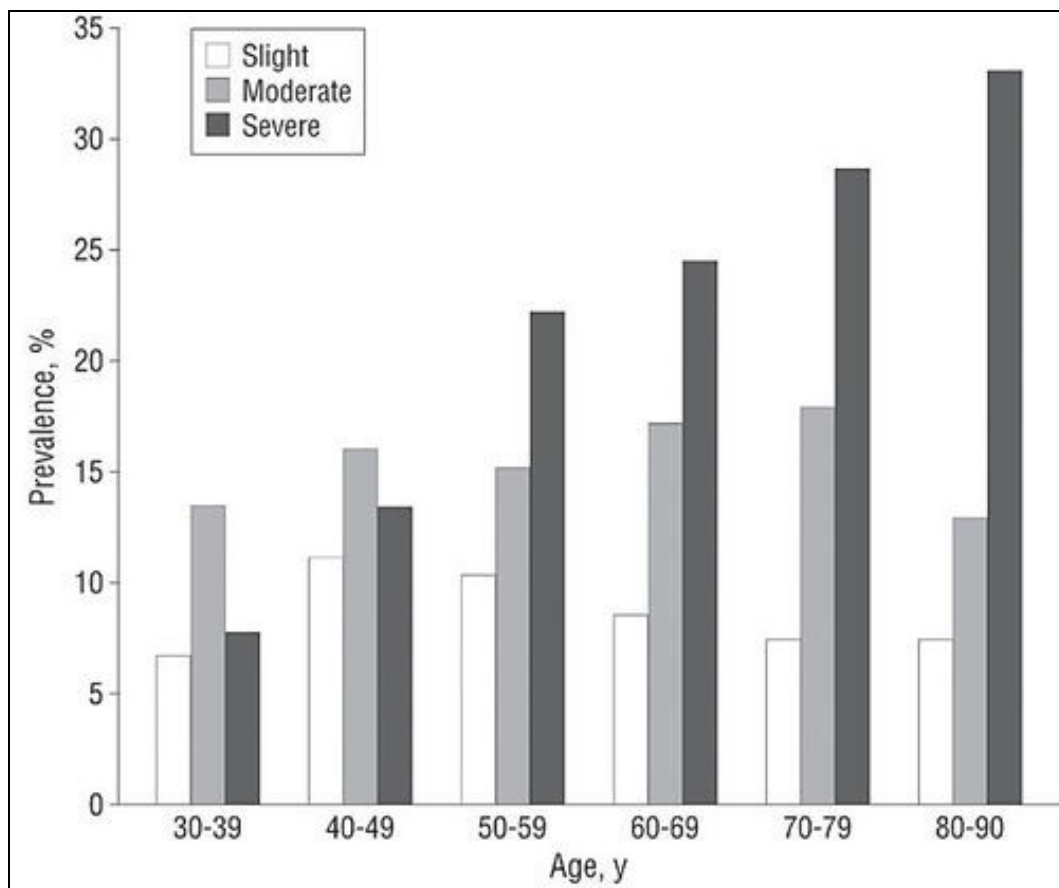


Source: Archives of internal medicine

### Market potential for ureteral stent in US (women)

Market potential for ureteral stents in US women, 2009					
1. Age groups	2. Female population (from US census data)	3. Prevalence rate in female (%)	4. Market potential (total prevalence) (2*3)	5. Catheterization rate (%)	6. Stent market based on catheterization rate (4*5)
30-39	20128402	28%	5635953	0.043%	2423
40-49	22074384	41%	9050497	0.123%	11132
50-59	20929761	48%	10046285	0.124%	12457
60-69	14605565	51%	7448838	0.160%	11918
70-79	9046207	55%	4975414	0.172%	8558
? 80	7216598	54%	3896963	0.044%	1715
<b>Total</b>			<b>41053950</b>		<b>48203</b>

- Prevalence rate in US women is growing at a CAGR of 1.19%
- Ureteral stent market is growing at a CAGR of 1.47%

### Ureteral stent market forecast in US (women)

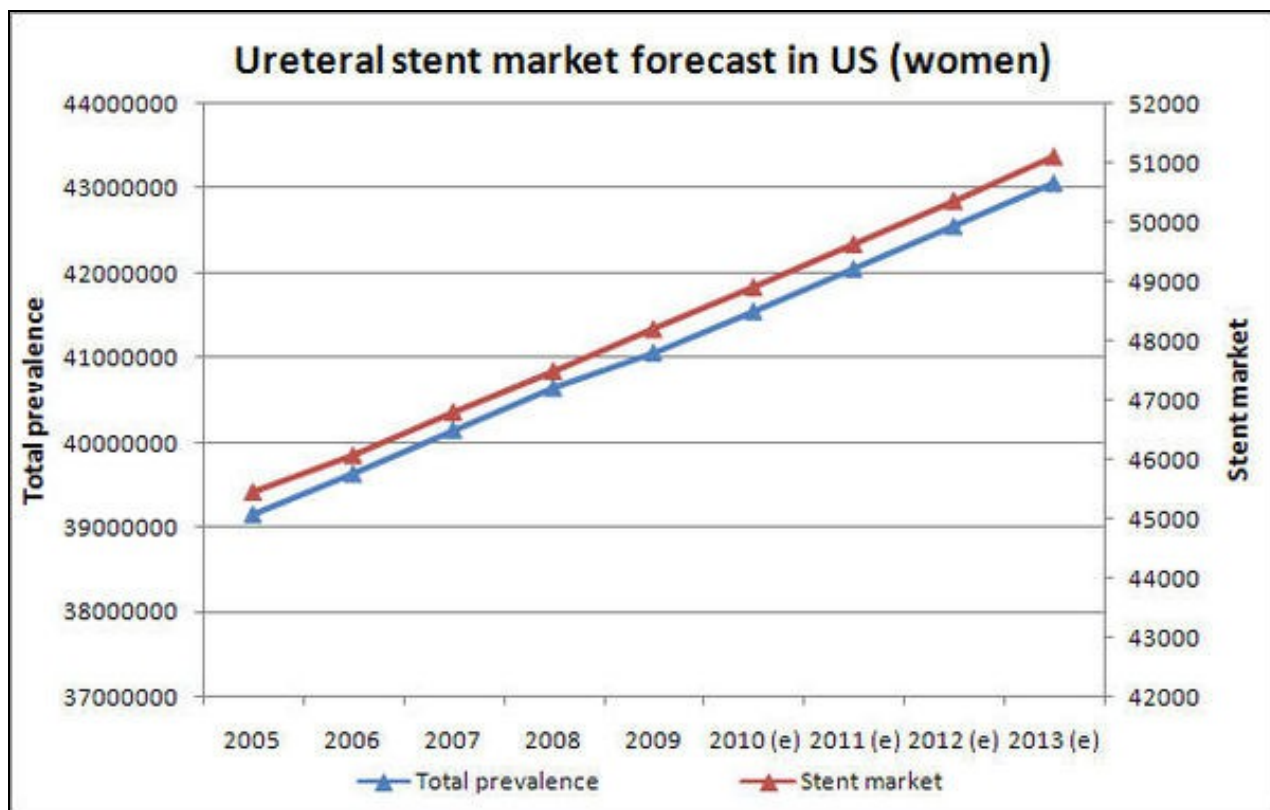

[Detailed model workbook](#)

## Ureteral stent companies

Various companies offering ureteral stents are:


- Boston Scientific Corporation
- Cook Group
- Applied Medicals
- Redi-Tech Medical Products

## Ureteral stents of various companies

### Boston Scientific


- Contour VL Variable Length Percuflex Stents

Inward spiral design of Nautilus Coil minimizes tissue contact for enhanced comfort. Dual variable length coil geometry permits balanced stent positioning to minimize the risk of migration. HydroPlus Coating provides unequalled surface lubricity which reduces friction to minimize risk of buckling during introduction and placement and to reduce risk of trauma and encrustation.

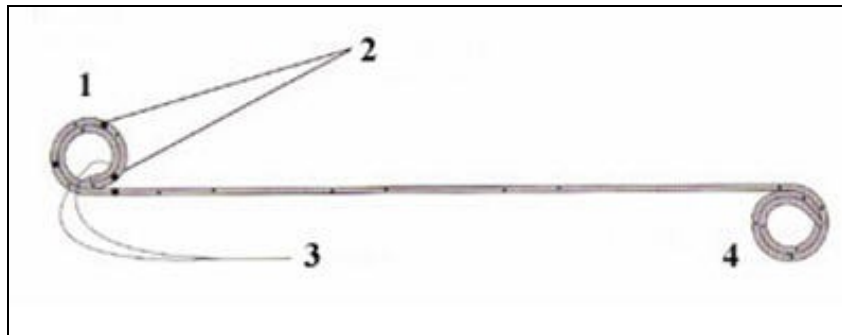

Source: [www.bostonscientific.com](http://www.bostonscientific.com)


- Percuflex Stents
  - ◊ High coil strength pigtail shape prevent stent migration
  - ◊ Thin wall promotes drainage and patency
  - ◊ Multiple, large side ports promote drainage
  - ◊ Attached suture for positioning and subsequent removal without the need for repeat cystoscopy

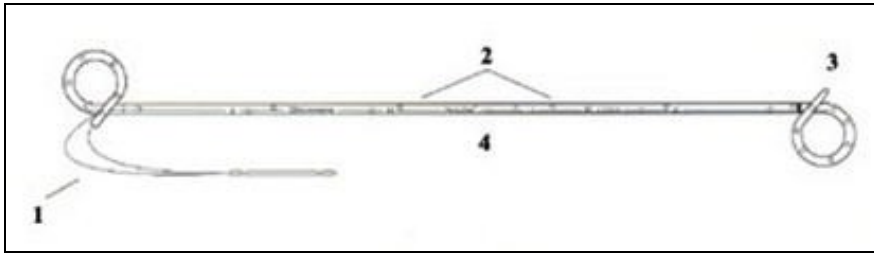

Source: [www.bostonscientific.com](http://www.bostonscientific.com)

- Polaris Ultra Ureteral Stent

- ◊ Ultra Ureteral Stent provides the ease of placement benefits of a firm stent graduated into a soft bladder coil
- ◊ Co-extrusion combines a firm durometer Percuflex Plus Material and a soft durometer Percuflex Material in the same stent
- ◊ Nautilus bladder coil potentially reduce bladder irritation, and a relaxed renal coil facilitate ease of removal

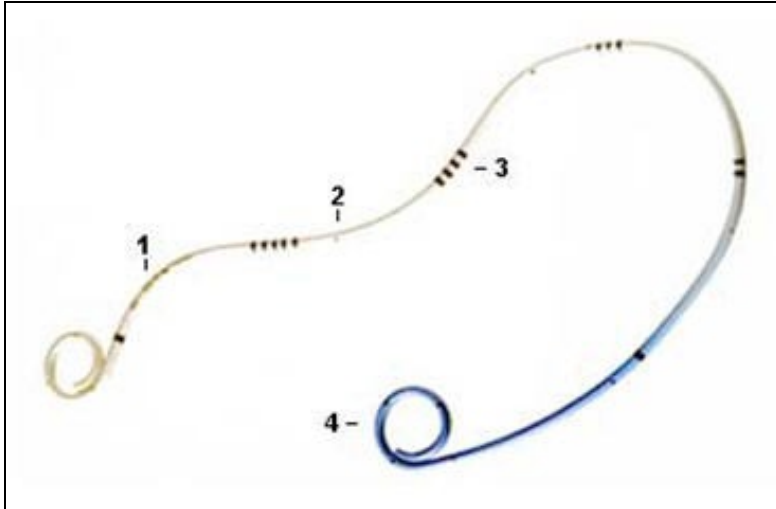

Source: [www.bostonscientific.com](http://www.bostonscientific.com)

Cook Group

- Bander Ureteral Diversion Stent Set


It is used for intraoperative placement to stent the ureter during ureteroileal conduit construction and continent urinary diversions. Set includes: 2 stents, 2 catheter retainers and wire guide.



Source: [www.cookmedical.com](http://www.cookmedical.com)

- C-Flex Double Pigtail Ureteral Stent Set

It is used for temporary internal drainage from the ureteropelvic junction to the bladder. Set includes stent, wire guide, stent positioner and catheter.

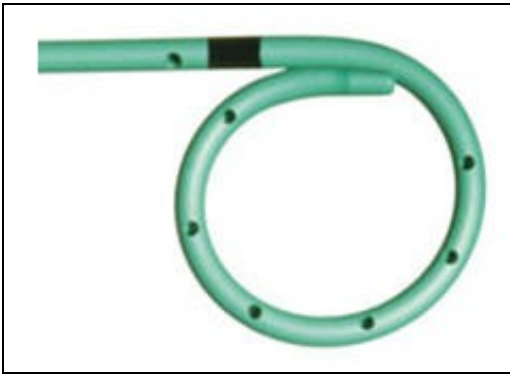

Source: [www.cookmedical.com](http://www.cookmedical.com)

- Towers Peripheral Ureteral Stent Set

It is also used for temporary internal drainage from the ureteropelvic junction to the bladder. The stent configuration allows peripheral as well as luminal drainage. Set Includes: Stent, Wire Guide, Catheter, and Stent Positioner.



Source: [www.cookmedical.com](http://www.cookmedical.com)

**Applied medicals**

- 7-10 endopyelotomy stent

It is used by urologists for endopyelotomy and endoureterotomy. The dual diameter promotes optimal healing while minimizing the discomfort often associated with larger diameter stents.




Source: [www.appliedmed.com](http://www.appliedmed.com)

**Redi-Tech Medical Products**

- Ureteral stents set

- ◊ Attached suture for positioning and subsequent removal without the need for repeat cystoscopy
- ◊ Multiple, large side ports promote drainage
- ◊ Radiopaque stent markings aid in placement and sizing

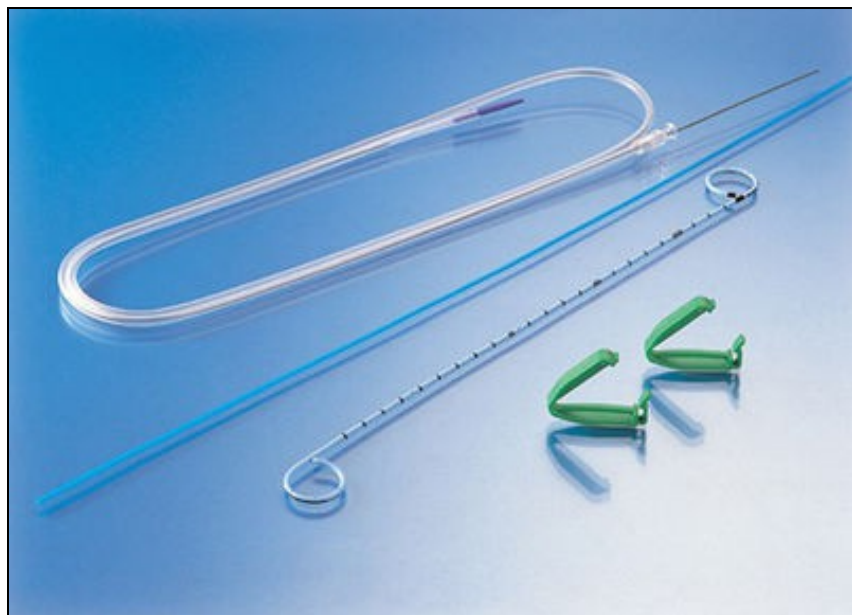

Source: [www.redi-tech.com](http://www.redi-tech.com)

**Like this report?**

**This is only a sample report with brief analysis  
Dolcera can provide a comprehensive report customized to your needs**

Buy the customized report from Dolcera		
Patent Analytics Services	Market Research Services	Purchase Patent Dashboard
Patent Landscape Services	Dolcera Processes	Industry Focus
Patent Search Services	Patent Alerting Services	Dolcera Tools

## Clinical Trials

### New trials

Title	Conditions	Intervention	Sponsors and Collaborators
Assessment of Drug-Eluting Ureteral Stent on Bacterial Adherence and Biofilm Formation	Renal Calculi, Ureteral Obstruction	Device: Ureteral Stent	Lawson Health Research Institute, Boston Scientific Corporation
Memokath® 044TW Stent for Treatment of Urethral Stricture	Urethral Stricture	Device: Memokath stenting	Engineers & Doctors Wallsten Medical Group
Study to Determine if There Are Specific Clinical Factors to Determine Stent Encrustation	Kidney Stones	N/A	University of California, Irvine
Ureteral Stent Length and Patient Symptoms	Kidney Stones	Device: Ureteral Stent	Emory University
Drainage of Malignant Extrinsic Ureteral Obstruction Using the Memokath Ureteral Stent	Ureteral Obstruction	Device: Memokath 051 Ureteral Stent	Mayo Clinic Engineers & Doctors Wallsten Medical Group
A Prospective Comparison Between Ureteral Stent and Nephrostomy Tube for an Urgent Drainage of Obstructed Kidney (JJVsPCN08)	Kidney Disease	Device: nephrostomy tube and ureteral stent	Rabin Medical Center

### Concluded trials

Title	Abstract	Enrollment	Disorder	Conclusion
Long-term outcome of permanent urethral stents in the treatment of detrusor-sphincter dyssynergia	To evaluate the long-term efficacy of a permanently implanted urethral stent in the treatment of spinally injured patients with detrusor-sphincter dyssynergia.	13	Detrusor-sphincter dyssynergia	Stenting is an effective alternative to sphincterotomy in the long-term, although secondary bladder neck obstruction is a frequent problem.
Nephrostomy Tube or 'JJ' Ureteric Stent in Ureteric Obstruction: Assessment of Patient Perspectives Using Quality-of-Life Survey and	Upper urinary tract obstruction is often relieved by either a percutaneous nephrostomy tube (PCN) or a ureteric stent. Both can cause considerable morbidity and	34	Upper urinary tract obstruction	Patients with 'JJ' stents have significantly more irritative urinary symptoms and a high chance of local discomfort than patients with nephrostomy tubes (PCN). However, based on the EuroQol analysis,

Utility Analysis	reduce patient's health-related quality of life (QoL). We have compared the QoL in these 2 groups.			there is no significant difference in the gross impact on the health-related QoL or the utility between these groups indicating no patient preference for either modality of treatment.
Impact of stents on urological complications and health care expenditure in renal transplant recipients: results of a prospective, randomized clinical trial.	A randomized, prospective trial to compare the incidence of early urological complications and health care expenditures in renal transplant recipients with or without ureteral stenting.	201	Renal transplant recipient	Using a ureteral stent at renal transplantation significantly decreases the early urinary complications of urine leakage and obstruction. However, there is a significant increase in urinary tract infections, primarily beyond 30 days after transplantation. Stent removal within 4 weeks of insertion appears advisable.

#### Adverse Events

S. No.	Brand Name	Adverse Event	Date FDA Received
1	<a href="#">Cook Urologicals Cook Urological Stent</a>	Stent broke into pieces while removing it from the patients body.	12/14/2005
2	<a href="#">Boston Scientific Boston Scientific Ureteral stent System</a>	Fractured stent seen under Fluroscopy	10/17/2007
3	<a href="#">Boston Scientific Boston Scientific Ureteral Stent System Kit 8 FR X 24 CM</a>	During insertion of ureteral stent, the stent broke into multiple parts which were retained in the patient.	10/14/2005
4	<a href="#">Boston Scientific Corp Boston Scientific 8 FR X 28 CM Ureteral Stent System Kit</a>	Breakage of the upper loop of the ureteral stent while trying to insert it.	1/5/2005
5	<a href="#">Boston Scientific Bostoon Scientific Micro Vasive Contour VL Ureteral Stent</a>	Broken stent observed during x-ray procedure.	12/12/2003

#### Review Articles Non Patent Analysis

#### Products

Boston Scientific Scimed, Inc.	Cook Urological Incorporated	OptiMed Global Care
Polaris? Ultra Ureteral Stent	Firlit-Kluge Urethral Stent	Opti-J Ureteral Stent System
Polaris? Loop Ureteral Stent	Koyle Diaper Stent	Ureteral Stent Sets, ureterorenoscope
Stretch? VL Variable Length Flexima® Stents	Silicone Universal Drainage Stent	Extra Strong Stent Sets (-Tumor)
	Tarkington Urethral Stent Set	Steerable Ureteral Stent Sets
	Zaontz Urethral Stent	Multilength
	Pediatric Urethral C-Stent	

#### Startup activity

- AbbeyMoor Medical Inc., a med-tech firm that's developed devices for treating urological disorders, has raised \$2.7 million in bridge financing.

## Phase 2: Deeper Dive


### Scenario

Client wishes to acquire a ureteral stent company.

### Deal analysis for a target company



**Key Deal Model Implications**

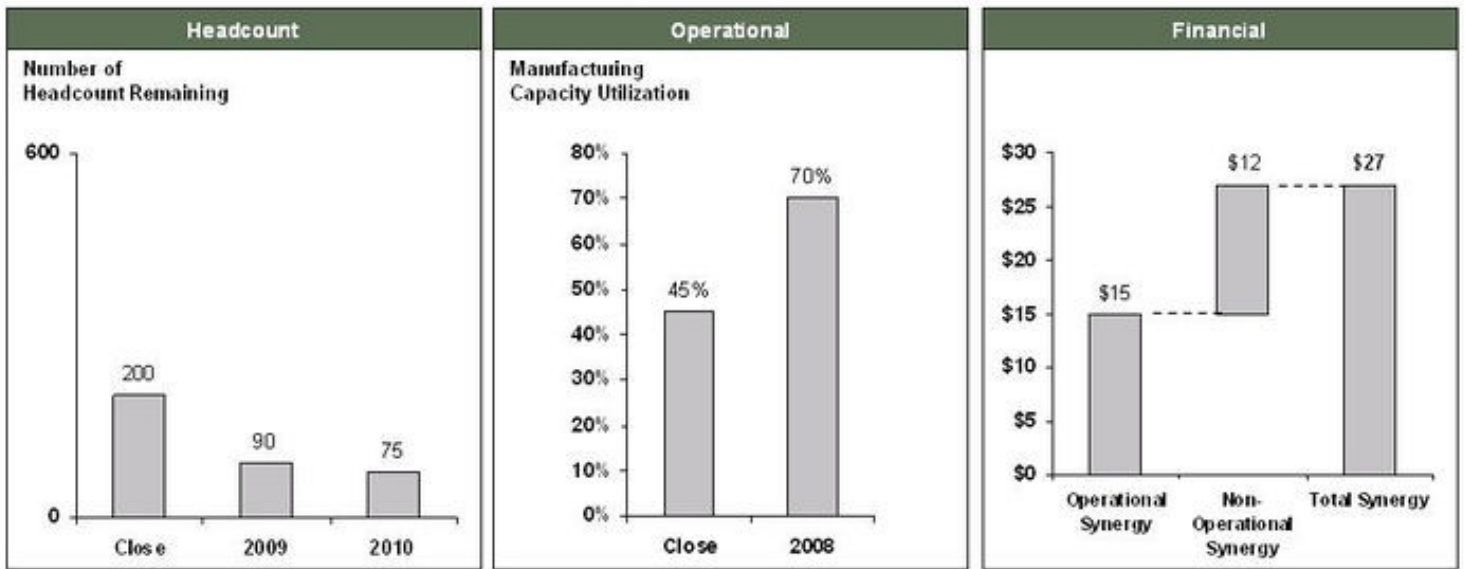

Deal implications

**Design History File Review: Review components**

Review	Verification	Tasks	Expertise
<b>Design Input</b>	Design input documents for sufficiency	<ol style="list-style-type: none"> <li>1. Access electronic data room.</li> <li>2. Check what documents are provided.</li> <li>3. Compare document list with standard client document list.</li> <li>4. Check whether each specified document has appropriate content.</li> </ol>	Quality systems
	Design input documents linked to the product performance specifications	<ol style="list-style-type: none"> <li>1. Compare product specifications to design inputs</li> <li>2. Check whether appropriate verifications and validations are performed</li> <li>3. Establish if all specifications are linked to design inputs</li> </ol>	
<b>Product Performance Specifications (PPS)</b>	Design inputs correlate adequately to the specifications; DV&V (design verification and validation) criteria are based on risk management documentation or if the criteria are based on sound statistical sampling plans	<ol style="list-style-type: none"> <li>1. Compare product specifications to design inputs</li> <li>2. Check whether appropriate verifications and validations are performed</li> <li>3. Establish if all specifications are linked to design inputs</li> </ol>	Quality systems, CAD
	Appropriate design verification and validations (DV&V) are performed	<ol style="list-style-type: none"> <li>1. Show DV&amp;V criteria are based on risk management requirements</li> </ol>	
	Product performance specifications correspond to appropriate design output documents	<ol style="list-style-type: none"> <li>1. Correlate design drawings with the specifications</li> <li>2. Check whether maximum dimensions, sizes etc. (with tolerances) are within the specified range</li> </ol>	
<b>Risk Management Documents</b>	Risk Analysis, Design Failure Modes and Effects Analysis (DFMEA), Process FMEA, other risk management documentation	<ol style="list-style-type: none"> <li>1. Check whether documentation is available</li> <li>2. Check whether it adheres to appropriate ISO 14971 standards</li> <li>3. Check whether it adheres to appropriate client standards</li> </ol>	Quality systems
	DFMEA links appropriately to the PPS	<ol style="list-style-type: none"> <li>1. Verify whether DFMEA and product specifications are correlated</li> </ol>	
	Appropriate DV&V reports and design output documents are referenced correctly as risk mitigation activities in the DFMEA	<ol style="list-style-type: none"> <li>1. Validate the process and correlate with design inputs</li> <li>2. Validate that sizes used are within range of risk mitigation criteria</li> </ol>	

	PFMEA links appropriately to the process validation protocol acceptance criteria; In-process inspection procedures and/or manufacturing procedures are recorded as appropriate risk mitigation activities in the PFMEA	1. Validate the process protocol 2. Validate the inspection procedures used	
Design Output Documents	Completeness of drawings	1. Check if the CAD diagrams overlay and "fit" perfectly 2. Check tolerance stackups	Quality systems, CAD
	Correlate First Article Inspection data to the dimensions on the drawings	1. Obtain First Article Inspection data 2. Check if this data correlates with the completeness of drawings	
Manufacturing Documents	Manufacturing procedures, component specifications, raw material specifications, incoming and in-process inspection procedures for completeness	1. Verify the Bill of Materials corresponds to raw materials and manufacturing procedures 2. Correlate incoming and in-process inspection procedures with the process specifications	Material science, manufacturing engineering, quality systems
	Linkage between component and raw material specifications and appropriate incoming inspection procedures	1. Identify any missing documentation for inspection procedures	
	Inspection procedures have adequate sampling plans based on PFMEA risk mitigation levels ? this includes packaging and labeling materials	1. Review supplier audit reports for compliance	
	Calibration records and preventive maintenance records; in-process / incoming inspection test methods and related test method validations	1. Check the entire equipment-related lifecycle 2. Check if machine operational qualification was performed 3. Check if the measurement equipment was validated	
Validation Report	DV&V reports, Shelf-life reports, Biocompatibility test reports, Sterilization reports, Packaging Validation reports, Process Validation Reports	1. Ensure all reports are available and linked together appropriately 2. Identify all inconsistencies across different reports	Quality systems
	Design test methods and related test method validations	1. Compare test methods used to those in client and ISO standards 2. Identify inconsistencies across test methods	

### Sample report

#### Performance/Functional Characteristics

Design Input				Design Output	Design Verification Report #	Status (P/F/R)	Design Validation Report #	Status (P/F/R)
User Needs	User Need Rationale	Engineering Specification	Engineering Specification Rationale					
Provide antimicrobial resistance for up to 2 weeks	Ureteral Stent User Survey (Document #XXXXX)	Stent must have chlorohexadine surface concentration of 10-20 mg/cm <sup>2</sup> for 3 weeks	Document #XXXXX	Test Document #XXXXX	Report 01-005-06-007	P	Report 01-005-06-007	P

### Potential DHF Review Outcomes

Based on a review of the above DHF documents a potential outcome for the uretral stent acquisition project could involve the following:

1. Better explanation of existing design input documents and also better linkage between the design inputs and product specifications.
2. Creation of some new test methods for design, incoming and in-process inspections and also include recommendations for the test method validations. Creation of any new DV&V data would be highly unlikely as it could potentially trigger a new submission or a note-to-file to the regulatory agencies.
3. Change in raw materials to better grade materials e.g. Switching resin to a USP Class VI biocompatible resin. This would eliminate some on-going testing but require additional upfront one time biocompatibility testing.
4. Updating drawings based on results from the FAI data.
5. Converting existing Company Y documents into Company X format and identifying potential gaps and streamlining linkage between raw material specifications and inspection procedures.
6. Identifying installation, operational and process qualification requirements with the assumption that no additional design verification and validation activities are required based on the fact that the device is currently approved for sale in the US and ROW.
7. Recommend activities necessary for completing packaging, labeling, ship testing and shelf-life testing. Stress should be on being able to leverage existing data for shelf-life without changing the regulatory status of the device.
8. Company X may want to perform additional biocompatibility testing to create an internal baseline and also update their biocompatibility files.
9. Help streamline suppliers for components when switching over from Company Y to Company X. Search for existing Company X suppliers that can supply off the shelf items that Company Y may be sourcing from other vendors / suppliers.
10. Identify process improvements that can be rolled into the manufacturing transfer without changing the design and impacting the existing regulatory status for the device e.g. instead of hand mixing pigment to resin use a pre-mixer to control quality of mixing and resulting extrusion

or perform the molding and over-molding steps in 1 machine instead of 2 separate molding machines.

## Phase 3: Post-acquisition integration

### Deadlines

**Goal:** Switch production transparently to new facilities transparently to the distribution system

Stage	Tasks	Milestone payment	Date
Design center integration plan	<ul style="list-style-type: none"> <li>Gap analysis completion (acquiree)</li> <li>Gap analysis completion (acquirer)</li> </ul>		September 15, 2008
Design to manufacturing transfer			December 15, 2008
Equipment transfer		Milestone I payment	Jan 7, 2009
Shut production at acquiree facility	Negotiation for contract extension	Milestone II payment	Feb 15, 2009
Start production in acquirer facility			Feb 7, 2009
Switch to new SKU			Feb 15, 2009
End development of new generation product/s in old facility			Feb 7, 2009
Restart development of new generation product/s post-acquisition		Final milestone payment	Mar 1, 2009

### Documents and Ownership

Document	Owner	Last update date
Product performance specifications	Paul Swain	07/27/2008 08:15:35 PST
Component specifications	Kevin Teller	06/12/2008 12:22:07 PST
Preclinical test results	Joanne Krannert	07/03/2008 14:17:00 PST
Clinical tests	Joanne Krannert	08/01/2008 08:00:55 PST

## Ureteral stents regulatory issues

The FDA classifies a ureteric stent as follows:

- TITLE 21 - FOOD AND DRUGS
- CHAPTER I - FOOD AND DRUG ADMINISTRATION DEPARTMENT OF HEALTH AND HUMAN SERVICES
- SUBCHAPTER H - MEDICAL DEVICES
- PART 876 - GASTROENTEROLOGY-UROLOGY DEVICES
- Subpart E - Surgical Devices
- Sec. 876.4620 - Ureteral stent.
- Classification - class II device [Code of Federal Regulations](#)

### Pre-Market Notification

Section 510(k) of the Food, Drug and Cosmetic Act requires device manufacturers who must register, to notify FDA of their intent to market a medical device at least 90 days in advance, also known as Premarket Notification. This premarket submission demonstrates to the FDA that the device to be marketed is at least as safe and effective, that is, *substantially equivalent*, to a legally marketed device. Parties required to submit a 510(k) to the FDA include domestic or foreign manufacturers introducing a device to the U.S. market, as well as specification developers and repackers/relabelers.




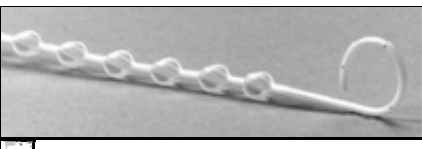
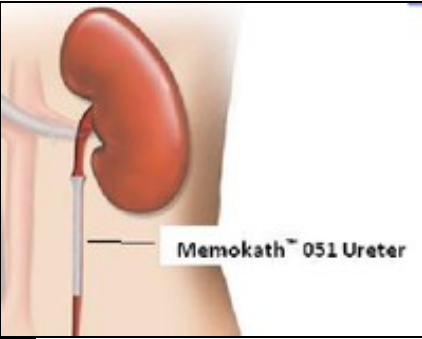
A 510(k) is required when:

- Introducing a device into commercial distribution (marketing) for the first time.
- Proposed different intended use for a device already in commercial distribution.
- Change or modification of a legally marketed device.

### 510(k) ?Substantial Equivalence? Decision Making Process

Some of the companies active in the field of ureteral stents have been represented in the table below. (This is not an exhaustive list and is just a sample)

Sr. No.	Company	Device	Approval	Date of Approval	Material	Technology	Indwelling time (days)	Image
---------	---------	--------	----------	------------------	----------	------------	------------------------	-------

1	Bard Urological	InLay Optima	FDA 510(k)	Dec 2004	Silicone	Double pigtail with monofilament suture loop	365	 InLay Optima
2	Boston Scientific	Polaris Loop	FDA 510(k)	Mar 2003	Dual Durometer Percuflex with HydroPlus Coating	Bladder loop design	365	 Polaris Loop
3	Cook Medical	Resonance	FDA 510(k)	May 2007	Metal	Temporary stenting	365	 Resonance
4	Fossa Medical	Stone Sweeper	FDA 510(k)	Aug 2002	Polyurethane	Spiral radially expanding stent	13	 Stone Sweeper
			CE Mark	Sep 2005				
5	Pnn Medical A/S	Memokath 051	CE Mark	1995	Nickel-titanium shape memory alloy	Double fluted ended spiral stent	240	 Memokath 051

## Timeline Sheet

### Ureteral Stent Timeline

#### • References

- 1- [http://en.wikipedia.org/wiki/Demographics\\_of\\_the\\_United\\_States](http://en.wikipedia.org/wiki/Demographics_of_the_United_States)
- 2- [http://www.wrongdiagnosis.com/c/catheter\\_infection/stats.htm?ktrack=kcplink](http://www.wrongdiagnosis.com/c/catheter_infection/stats.htm?ktrack=kcplink)
- 3- <http://www.fda.gov/downloads/ScienceResearch/SpecialTopics/WomensHealthResearch/UCM247851.pdf>
- 4- <http://www.managementparadise.com/forums/principles-management-p-o-m/208329-swot-analysis-boston-scientific-corporation.html>

## Like this report?

**This is only a sample report with brief analysis  
Dolcera can provide a comprehensive report customized to your needs**

Buy the customized report from Dolcera		
Patent Analytics Services	Market Research Services	Purchase Patent Dashboard
Patent Landscape Services	Dolcera Processes	Industry Focus
Patent Search Services	Patent Alerting Services	Dolcera Tools

## Contact Dolcera

**Samir Raiyani**

**Email:** [info@dolcera.com](mailto:info@dolcera.com)

**Phone:** +1-650-269-7952

## Addition





### Patent Categorization: Interactive mind map linked to Dolcera Dashboard

- To access the Dashboard you have to signup. You can do so by clicking [here](#)
- Use the mouse(click and drag/scroll up or down/click on nodes) to explore nodes in the detailed taxonomy
- Click on the red arrow adjacent to the node name to view the content for that particular node in the dashboard
- Click on the "+" sign to zoom the mindmap and "-" sign to shrink the mindmap

### Product-Patent-Clinical Trials Mapping

- To access the Dashboard you have to signup. You can do so by clicking [here](#)
- Use the mouse(click and drag/scroll up or down/click on nodes) to explore nodes in the detailed taxonomy
- Click on the red arrow adjacent to the node name to view the content for that particular node in the dashboard
- Click on the "+" sign to zoom the mindmap and "-" sign to shrink the mindmap


### Product to Patent Mapping

S. No	Company	Product	Patent no.	Date of Publication	Assignee/Applicant	Title	Patent Figure	co
1	<a href="#">Boston Scientific</a>	<a href="#">Percuflex® Ureteral Stent</a>	<a href="#">US5401257</a>	3/28/1995	Boston Scientific	Ureteral Stents, Drainage Tubes and the like		M pr en st m th po su E
2	<a href="#">Boston Scientific</a>	<a href="#">Percuflex® Plus Ureteral Stent</a>	<a href="#">US6719804</a>	4/13/2004	Boston Scientific	Medical Stent and Related Methods		M : co m fir va se a
3	<a href="#">Boston Scientific</a>	<a href="#">Polaris Ultra stent</a>	<a href="#">US6719804</a>	4/13/2004	Boston Scientific	Medical Stent and Related Methods		M : co m fir va se a
4	<a href="#">Boston Scientific</a>	<a href="#">Polaris™ Loop Ureteral Stent</a>	<a href="#">US6991614</a>	1/31/2006	Boston Scientific	Ureteral Stent for Improved Patient Comfort		M du of




Screenshot for the product to patent mapping(Bard and Boston)

- Click [here](#) to download the excel file.
- Mapped Patent vs Not Mapped Patents

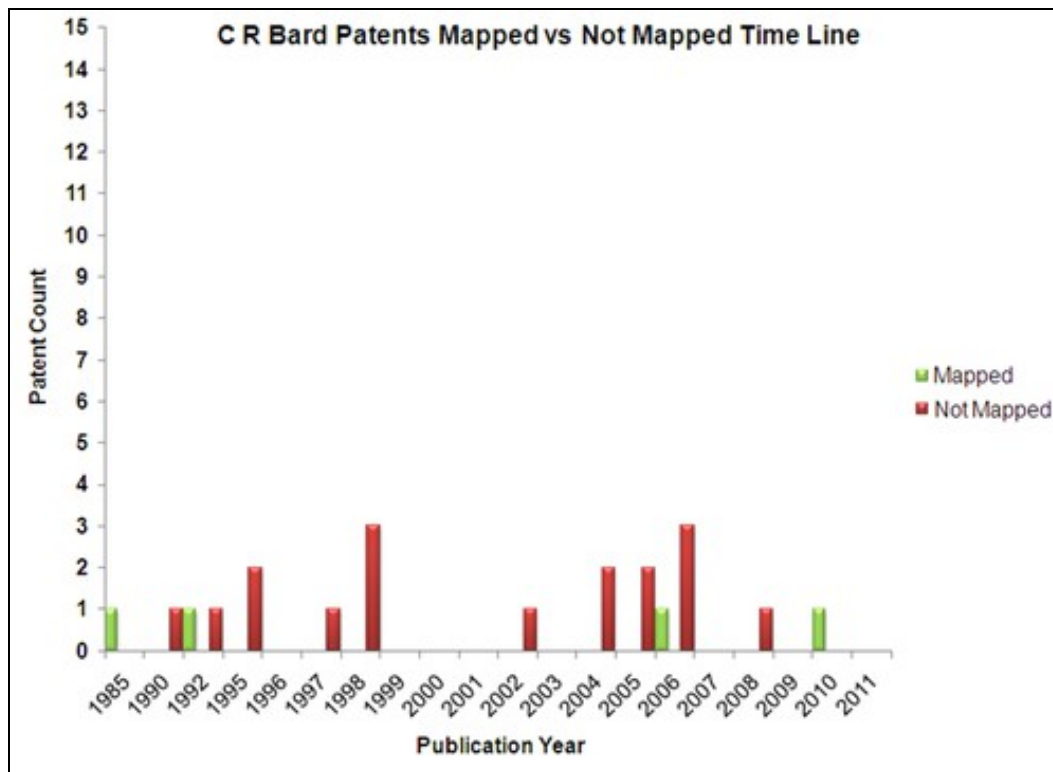

C R Bard

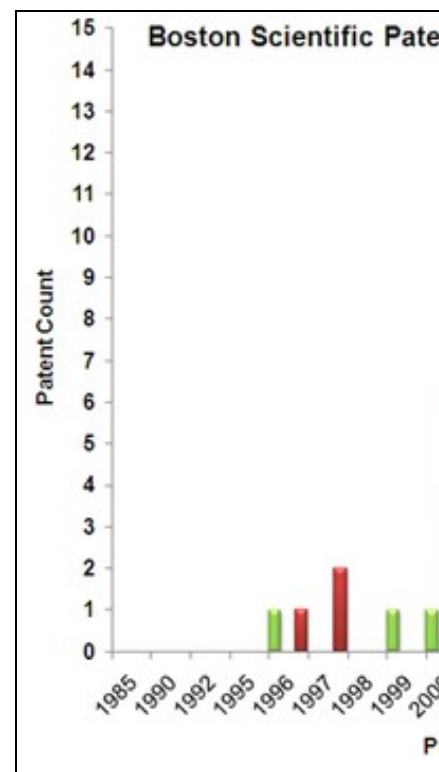

Boston Scientific



## Dolcera Dashboard

The screenshot shows the Dolcera Dashboard interface. On the left is a 'Data Filters' sidebar with a tree view of patent categories like 'Doubly fed induction generator', 'Parts (83)', 'Stator (69)', 'Rotor (73)', etc. The main area is titled 'Information' and contains a table of patent publications. Annotations include:

- Multi Level Classification**: A callout pointing to the 'Data Filters' sidebar.
- Assignees**: A callout pointing to a list of company names in the sidebar.
- Document PDF**: A callout pointing to a document icon in the patent table.
- Claims:**: A callout pointing to the 'Claims' section of a selected patent entry.
- Abstract:**: A callout pointing to the 'Abstract' section of a selected patent entry.
- Patent Charts**, **Patents**, and **Articles**: Callouts pointing to tabs at the top of the main area.
- Different Views**: A callout pointing to a view icon at the top.
- Export Selected Data**: A callout pointing to an export button at the top right.

Publication	Title	Assignee	Pub	Ap
US6448735B1	Controller for a wound rotor slip ring induction machine	Abb Research	2002	20
US20050189896A1	Method for controlling doubly-fed machine	Abb Research	2005	20
WO2007027141A1	Wind mill power flow control with dump load and power converter	Abb Research	2007	20
US20070114978A1	System for transmission of electric power	Abb Research	2007	20
US20090273187A1	Control method	Abb Research	2009	20
US20100085783A1	Method and system to influence the power generation of an adjustable speed generator	Abb Research	2010	20
US20090273187A1	Systems and methods for synchronous speed avoidance in doubly-fed induction generators	Acciona Windp	2009	20
US20100002475A1	Dc voltage regulator	Acciona Windp	2009	20
WO2010002402A1	Wind turbine control system and method	Acciona Windp	2010	20
GB2410386A	Low voltage ride through	American Supe	2010	20
GB2411252A	Low voltage ride through	American Supe	2010	20
GB2420456A	Generator control having grid imbalance detector	Areva T & D Uk	2005	20
US20100013343A1	Constant frequency and locked phase generator adaptable to variable torque	Areva T & D Uk	2005	20
US20030052643A1	Brushless doubly-fed induction machines employing dual cage rotors	Beijing Institut	2010	20
US20060192390A1	Method for operating a wind energy plant	Dewind Gmbh	2005	20
US20090021013A1	Brushless doubly-fed induction machine control	Dual Stator Te	2001	20
US20090302608A1	Control and protection of a doubly-fed induction generator system	Dual Stator Te	2003	20
US20090302608A1	Wind power system and method of operating it	Gamesa Innov	2006	20
US20090302608A1	Wind power installation and method of modifying the blade pitch in a wind power installation	Gamesa Innov	2009	20

### Dashboard Link

Ureteral Stent - Dashboard

- Flash Player is essential to view the Dolcera Dashboard
- To access the Dashboard you have to signup. You can do so by clicking [here](#)

## Key Artifacts

- Investment Heat Map
- Revenue Heat Map
- Patent Heat Map
- Dynamic Patent Dashboard
- Stent Landscape-Flash
- Company Profile - Flash

## Removed Sections

Removed Sections